EDGE^M TMS

Transcranial Magnetic Stimulation Chair

Innovative Medical Solutions

EDGE^M TMS

Transcranial Magnetic Stimulation Chair

EDGE^M TMS

A transformative and state of the art model for modern TMS (Transcranial Magnetic Stimulation) Therapy.

EDGE^M TMS is BTC ultimate concept Specialty Treatment Chair developed and manufactured exclusively for the most modern and efficient TMS Therapy both Clinical and Research.

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Registered Community Design

EDGE^M TMS

Thanks to an innovative engineering and the highest constructive approach EDGE^M TMS sets new and higher standards in the industry; the Memory Foam and independent padded sections are electrically adjusted guaranteeing the correct positioning and the best comfort for the patient during the treatment session, while the Foot Control and the Articulated Headrest allow for perfect positioning of the area during the procedure.

The unique design of BTC Articulated Headrest grants extreme versatility and precision, the side bars are adjustable (up/down and left/right) and removable, padded for better comfort and manufactured in Aluminium (including the screws) with a marking hole in the terminal block (available version without marking block). The back bar allows for perfect head stability and it allows for up/down and forward/backward adjustment thanks to an innovative articulation.

All surfaces are easy to disinfect and the quality of all materials is guaranteed for 2 (Electronic parts) and 5 (Mechanical parts) years.

EDGE^M TMS is paramount in Clinical and Research activity in the modern Neurostimulation and Therapy Centres.

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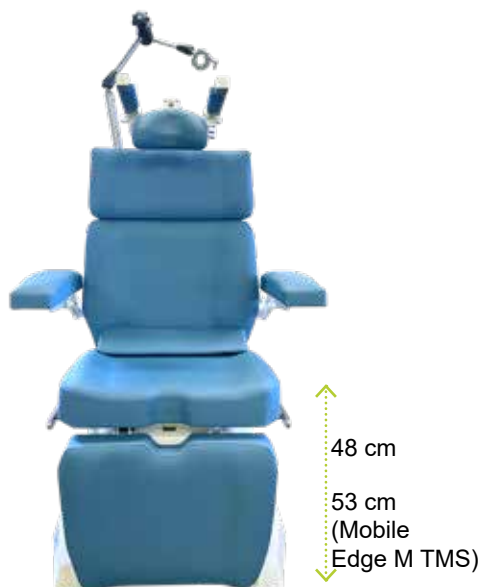
EDGE^M TMS

EDGE^M TMS chair reaches a minimum height of 48 cm allowing patients with different medical conditions to easily access the chair from the front or the sides without any barrier.

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EDGE^M TMS is equipped with medical grade castors, two of them are pivoting while 2 fixed, and central brake system at rear, in order to provide maximum maneuverability.

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360°

EDGE^M TMS can automatically reach the Operating Position thanks to the memorizable positions function activated by the foot control.

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EDGE^M TMS chair can reach a perfect flat position at any height, for the maximum comfort of the Operator while performing a clinical examination or professional procedure. Maximum height 88 cm.

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Zero position of the chair.

EDGE^MTMS chair automatically reaches "0" position to facilitate the access and the exit of the Patient.

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EDGE^MTMS allows the operator to perform "Synchro" movements of the Chair enables the combined and synchronized movement of the back-rest and leg-rest.

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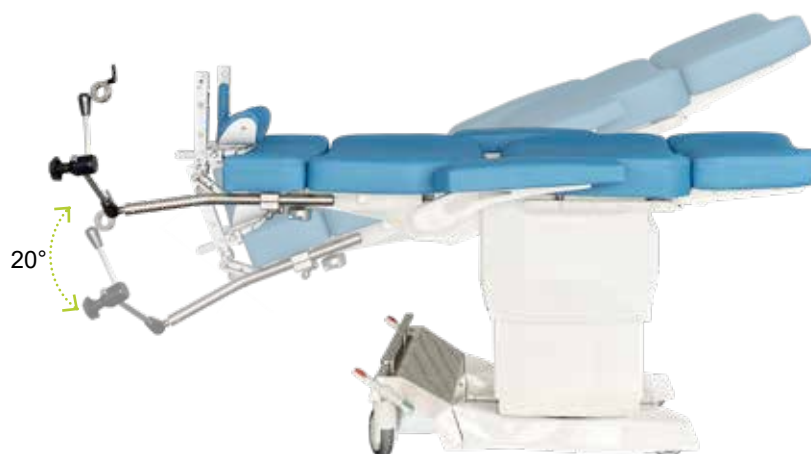
EDGE^MTMS can rely on a lifting capacity of 325 kg, one of the highest in its category.

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EDGE^MTMS reaches Trendelenburg position with 20° inclination.

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descrizione

9110666

ARTICULATED HEADREST WITH MARKING BLOCK

codice

descrizione

9110667

ARTICULATED HEADREST WITHOUT MARKING BLOCK

The side bars are equipped with a hole allowing the marking for the correct positioning of the Patient. Bars are also available without marker hole.

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The side bars are adjustable (up/down and left/right) and removable, padded for better comfort and manufactured in Aluminium (including the screws). The back bar allows for perfect head stability and it allows for up/down and forward/backward adjustment thanks to an innovative articulation.

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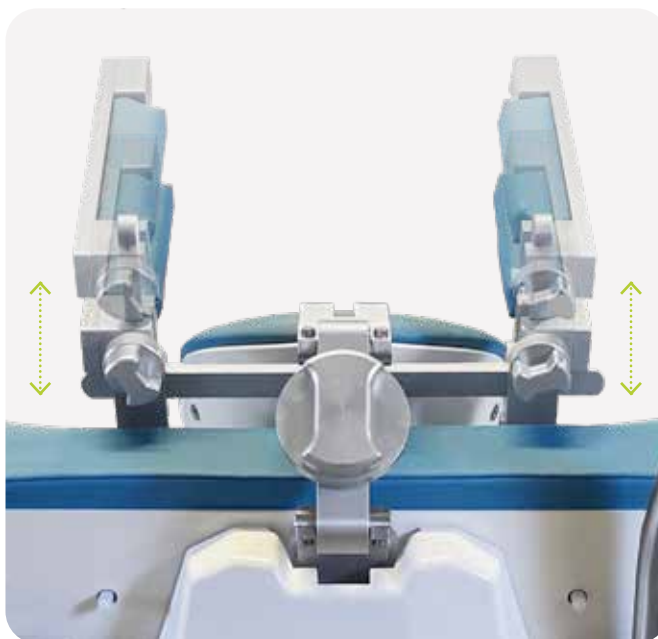
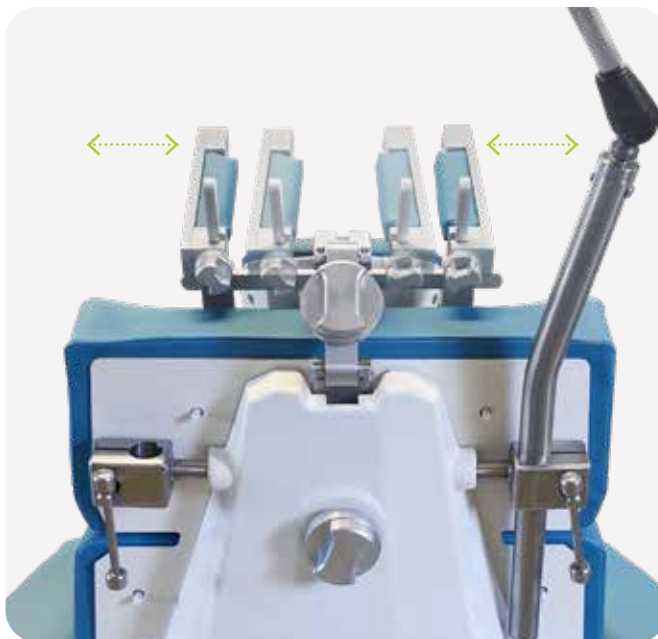
9110668

**COIL
SUPPORT**

The Coil support is equipped with an adjustable arm. The attachment is suitable with diameters between 28 and 34 mm. The support rod can be fixed both on the right side and on the left side.

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OPZIONALE

Arm-rest is shaped and padded according to the requirements of different medical sectors and can all be completely removed when necessary.

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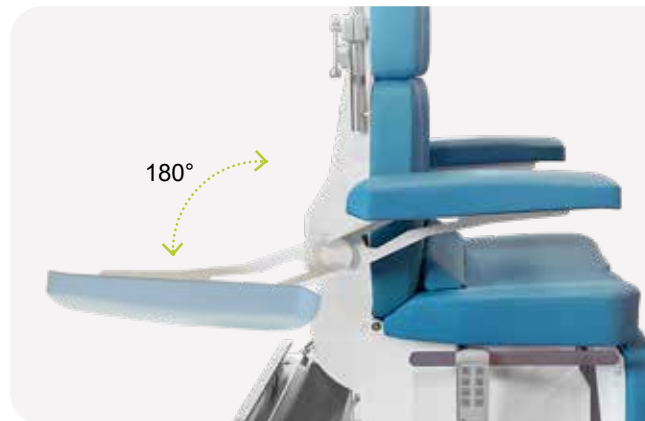
Both padded arm-rest are independently collapsable, removable and mechanically combined with the seat/back-rest movements ensuring comfort for the patient in any position.

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Padded and shaped leg-rest.
PROVIDED IN THE STANDARD VERSION

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Leg-rest operated with a multifunction foot control by with a linear actuator.

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2 Stainless steel side rails mounted on the seat.
PROVIDED IN THE STANDARD VERSION

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EDGE^M TMS is provided in the standard version with the structure in White color or Silver at choice.

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The upholstery comprises 4 separate padded sections fully detachable for easy cleaning and disinfection, made of high density undeformable foam and the cover is a bi-elastic, abrasion-resistant, flame retardant and antimicrobial medical grade Skai available in over 80 different colours and textures.

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EDGE^M TMS

EDGE^M TMS is operated by a wireless foot control with modern design and extremely intuitive. Easy to operate and clean, the foot control provides the Professional with an ergonomic and efficient working tool, useful to position the Patient comfortably and memorize different working or safety positions.

The foot control has 8 buttons to operate the 4 actuator and move the Chair sections; it is also provided with 2 memorizable positions (M1,M2), the easy Patient Access-Exit position (M0), and the Sync button that enables leg-rest and back-rest simultaneous movement.

The foot control is also equipped with two buttons for the emergency positions of the Chair: the flat position and the Trendelenburg position. The pedal innovative visual and acoustic signals communicate all safety, memorized movements and charging activity of the chair and it's foot control.

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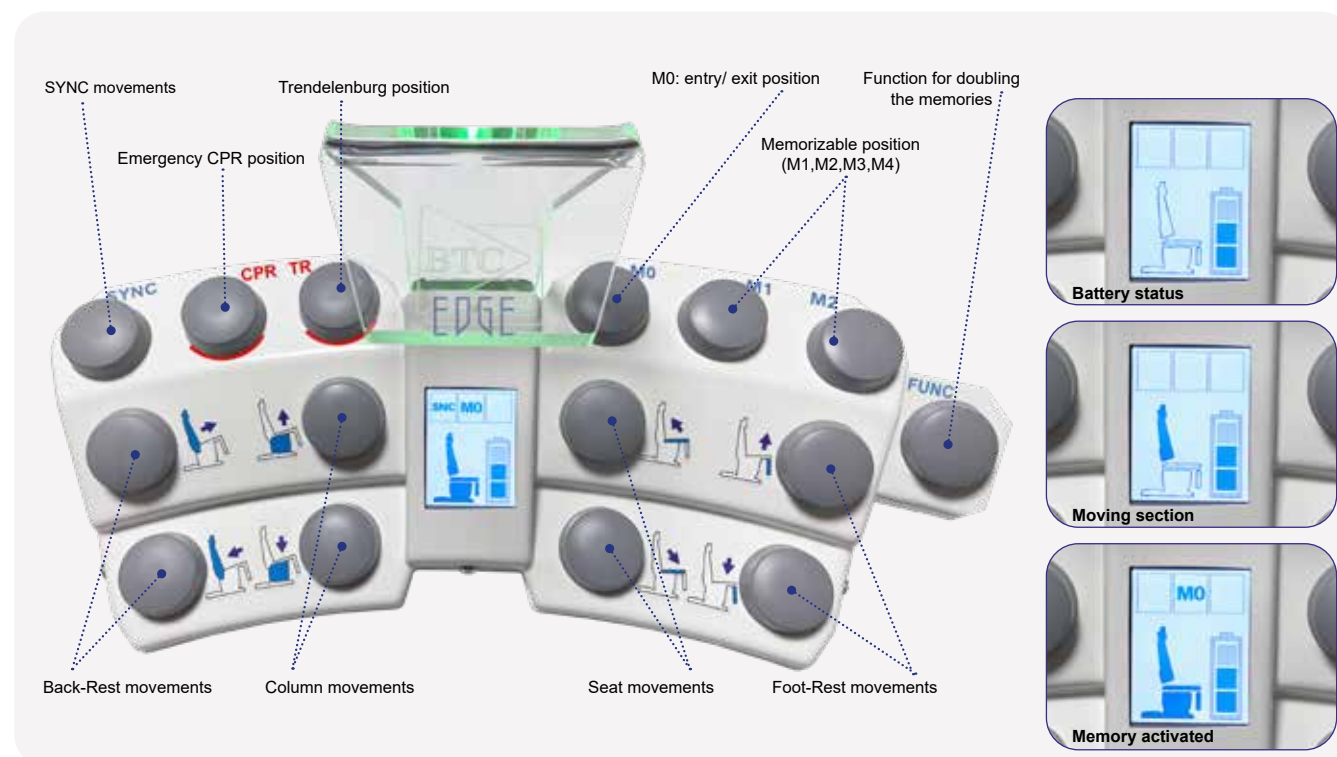
OPZIONALE

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The Hand Control enables the required positioning of the Patient, by activating the linear actuators that control the position of each independent section (lifting column, seat, back-rest and foot-rest).

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Body Velcro tape for Patient stabilization and containment.

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Belt with buckle for stabilization and containment of the Patient.

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Battery kit including two batteries (capacity 5Ah), with capacity indicator and recharger, allowing to operate the Chair without electrical supply connection. Replacing the battery with the spare one provided is very simple.

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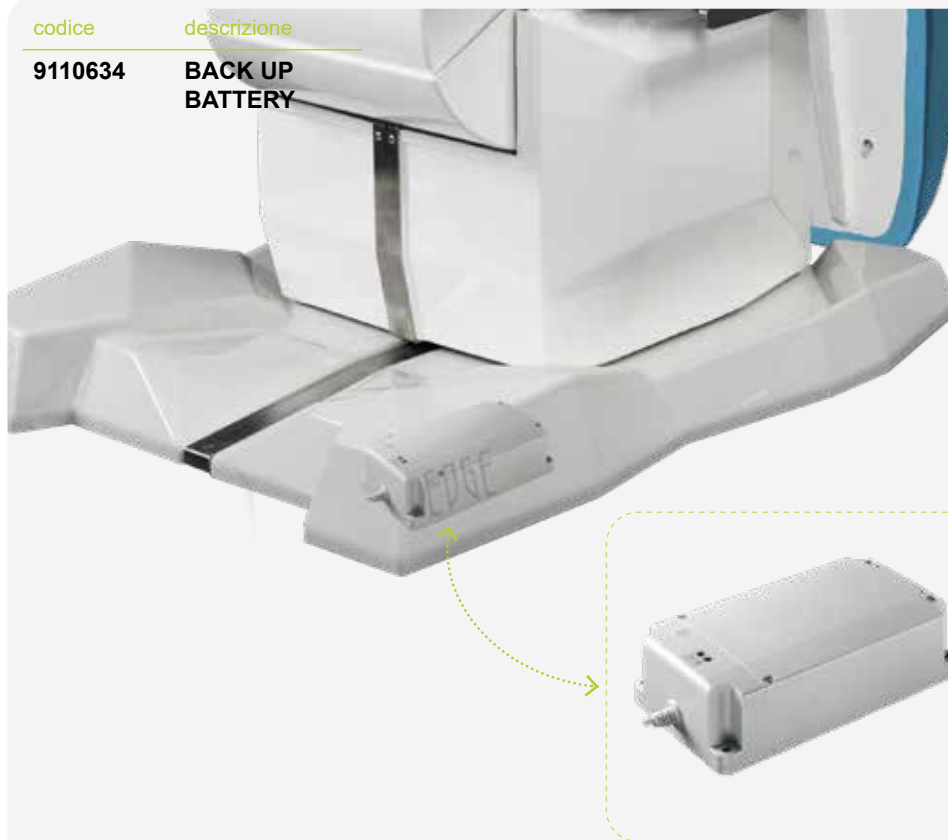
Back up battery with 2Ah capacity allowing to operate the Chair in the event of a temporary power outage.

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Foot control support for back rail.
OPTIONAL ON MOBILE EDGE M TMS

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codice	descrizione
9110631	RACK FOR MEDICAL ACCESSORIES (white structure)
9110631S	RACK FOR MEDICAL ACCESSORIES (silver structure)

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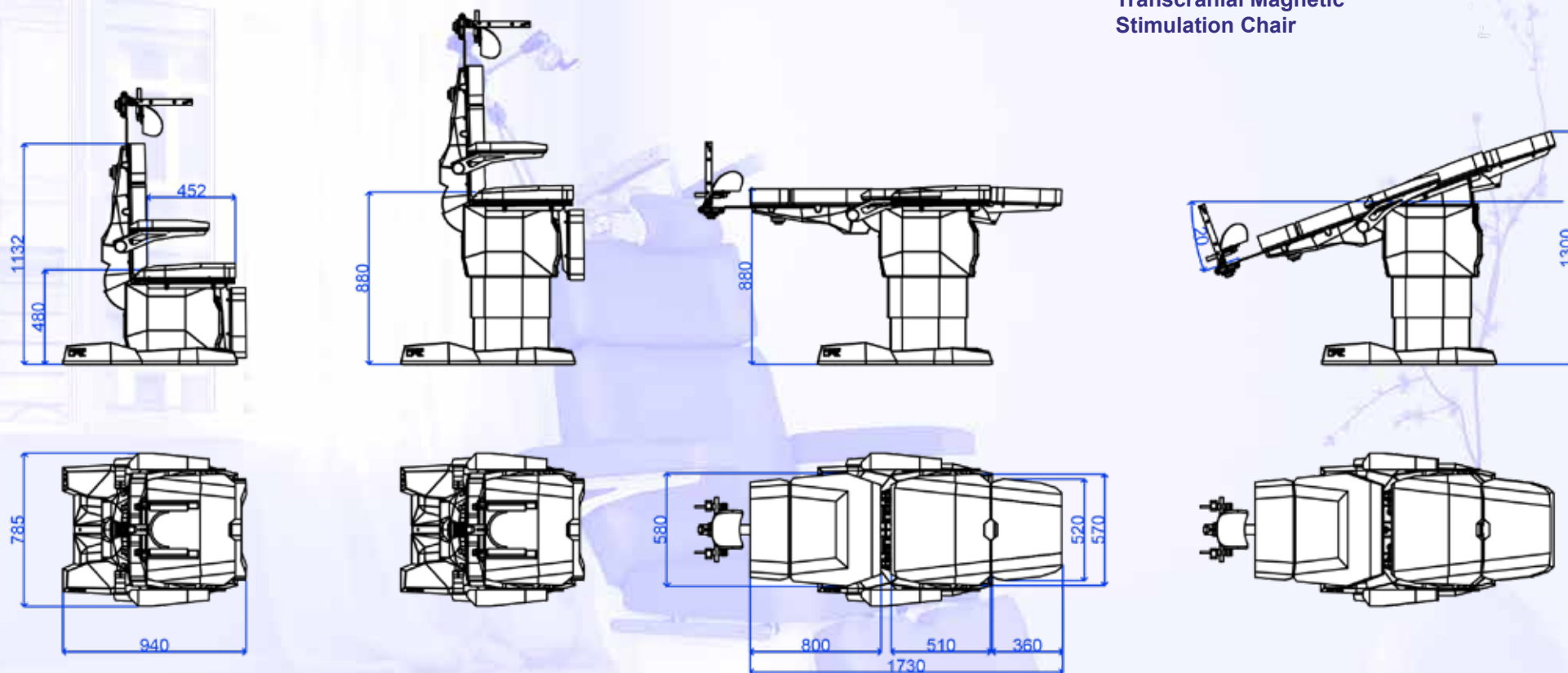
EDGE^M TMS

Transcranial Magnetic
Stimulation Chair

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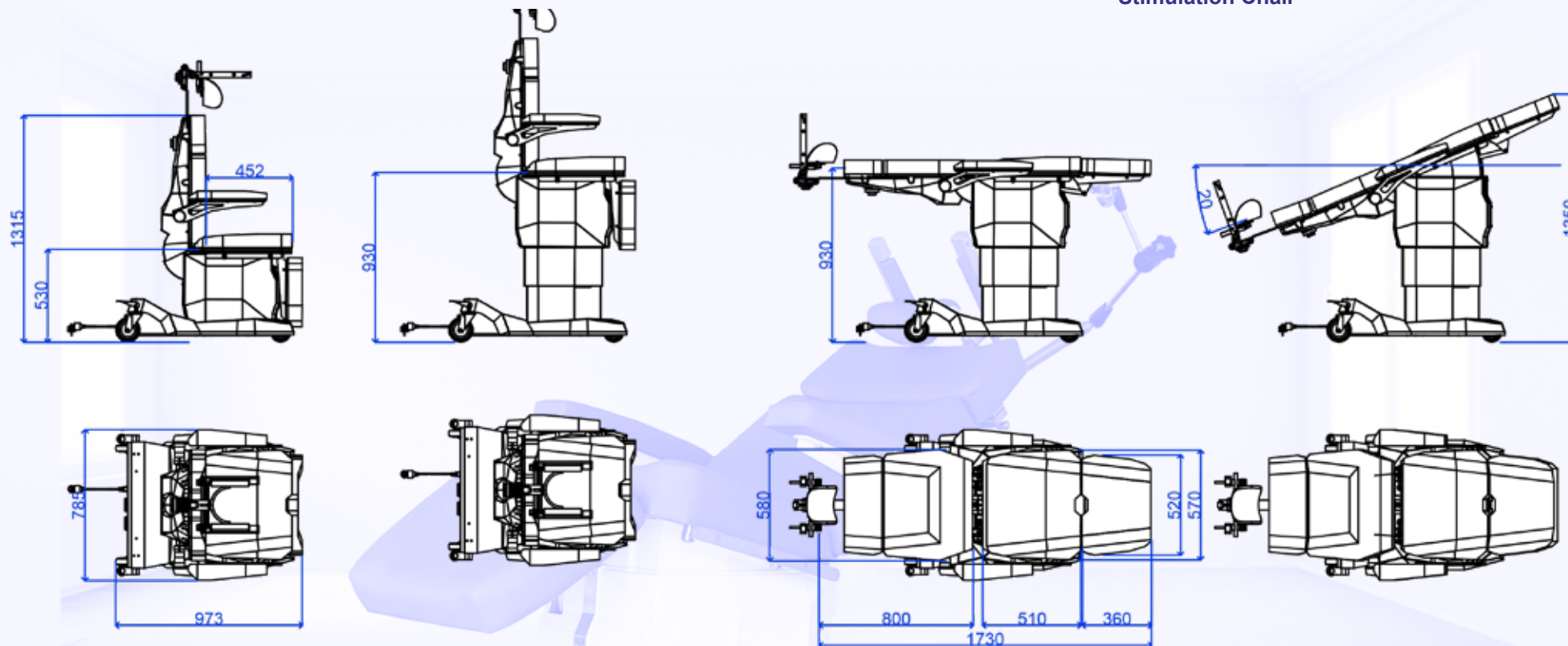
The seat naturally adapts to the user's body by increasing the contact area and decreasing the pressure on it. Circular seat with visco-elastic foam in different densities and memory effect. Multifunction armrest/back rest in high resistant polyurethane. Base in steel and polished chrome. 5 self braking wheels. Adjustable height by gas piston. High quality skay cover.

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Doctor stool with anatomical seat. Adjustable height of the backrest. Inclination angle of the backrest is adjustable in a fixed position or can dynamically support the back muscles. Control lever for all adjustments. Base in chromed steel. Seat height adjustable by gas spring. High quality skay coating.

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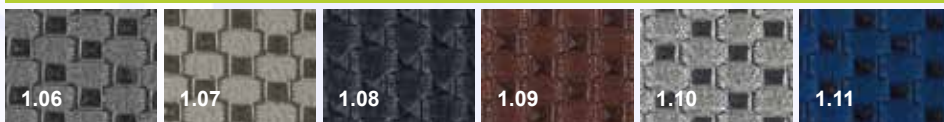
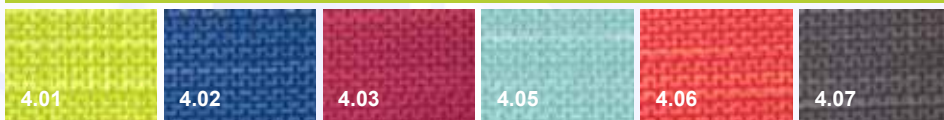
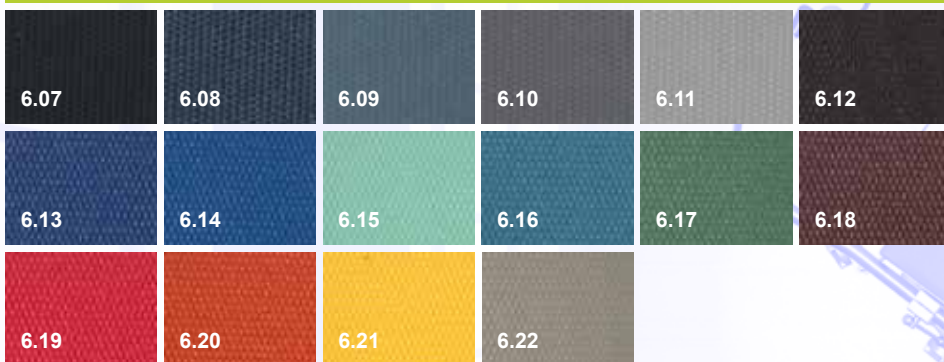
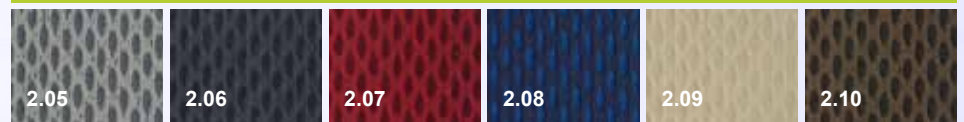
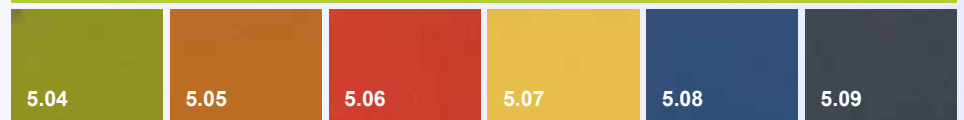
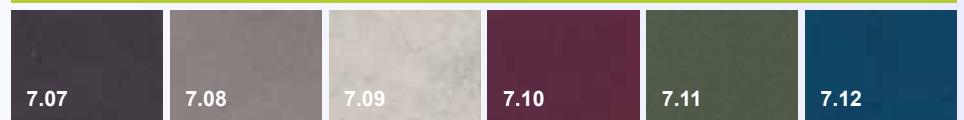
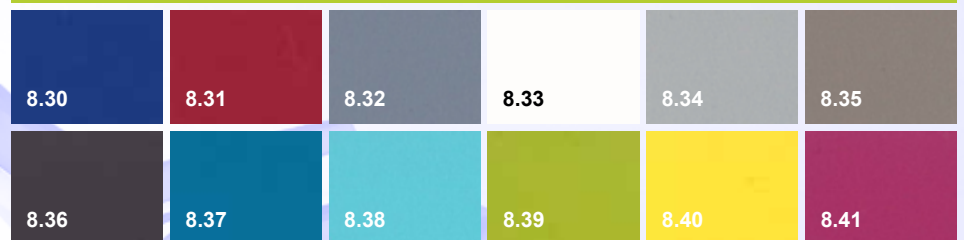
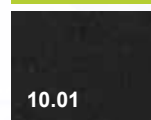
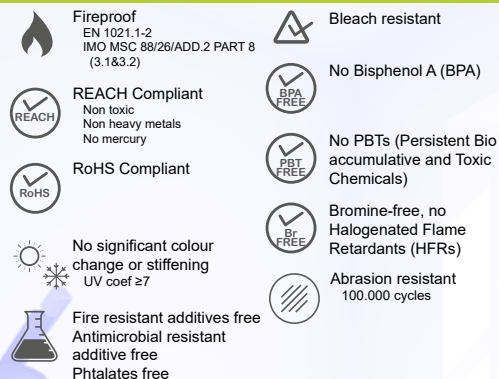
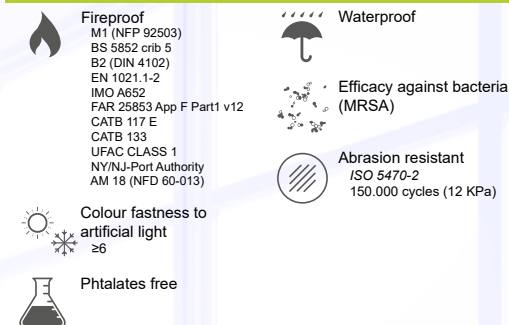
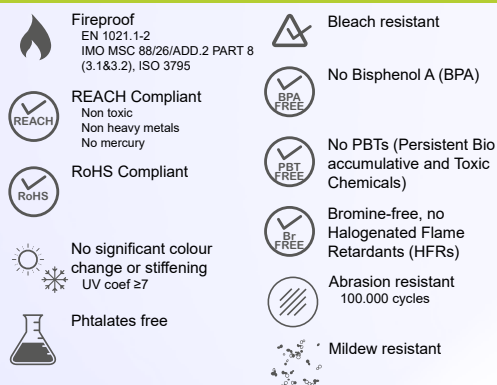
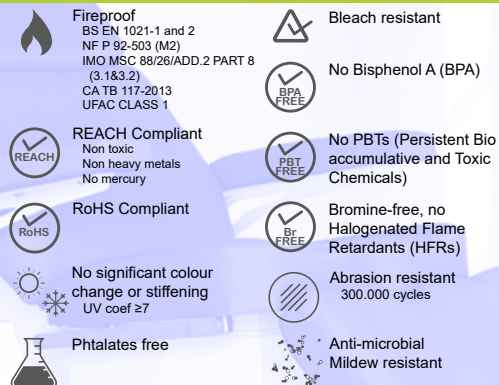
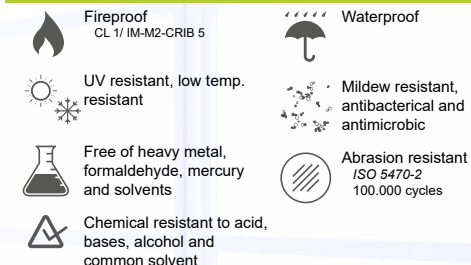
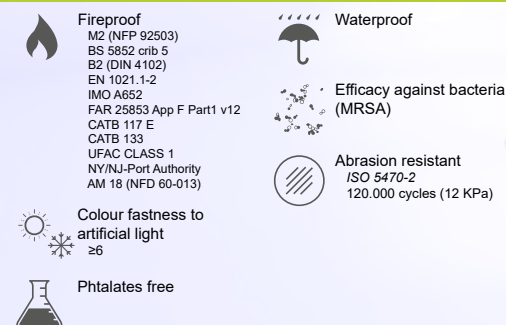
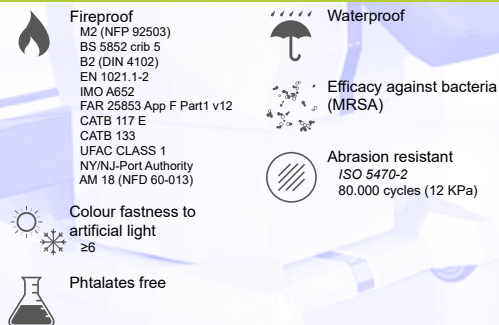
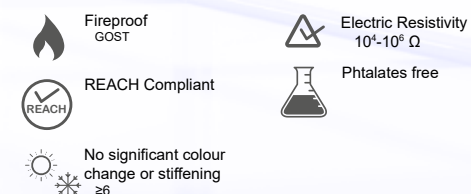
BTC Premium

BTC Antistatic

Electrically conductive (antistatic) synthetic leather

www.btc-med.it

btcmmed@btc-med.it

[btcmedicalequipment](https://www.facebook.com/btcmedicalequipment)

[btcmedicalequipment](https://www.instagram.com/btcmedicalequipment)

*"All BTC Medical Equipment is designed ,
developed and manufactured in Italy"*

USA Exclusive Distributor

www.jalimedical.com

inquiries@jalimedical.com

(781) 376 0008

BTC Srl

via Colombarotto, 34
40026 Imola (BO)
Tel. +39 0542 643664